

VAPORIZER PIN INDEX FILL VALVE UPGRADE

INSTALLATION PROCEDURE

1. Drain the vaporizer into a suitable container. Refer to the appropriate Narkomed Operator's Manual for draining instructions; **Observe all listed safety precautions when handling anesthetic agent.**

Refer to Figure 1 for the following steps:


2. Remove the two machine screws Ⓐ and remove valve retainer clip Ⓑ.
3. Unscrew existing valve Ⓒ and remove it from the block. Discard parts Ⓐ, Ⓑ and Ⓒ.
4. Verify that there are two seals in place on new valve assembly Ⓕ.

CAUTION

Ensure that the valve stem in new assembly Ⓕ is turned completely counter-clockwise with regard to the threaded hex gland.

5. Screw valve assembly Ⓕ into the vaporizer block, making sure that the valve stem position is maintained as explained in the above caution.
6. Carefully tighten the hex gland in the vaporizer block with a 9/16 in. open end wrench.
7. Close the valve by turning its stem clockwise.
8. Remove machine screw Ⓓ, install chain Ⓔ on the screw and reinstall Ⓓ into the vaporizer block. (The chain was previously attached by one of the screws on top of the block.)
9. Loosen screw Ⓖ securing the other end of the chain to the filler port plug. Turn the chain coupling to point downward, and then tighten the screw.
10. Perform a complete PMS procedure on the machine.

INSTALLATION PROCEDURE (continued)

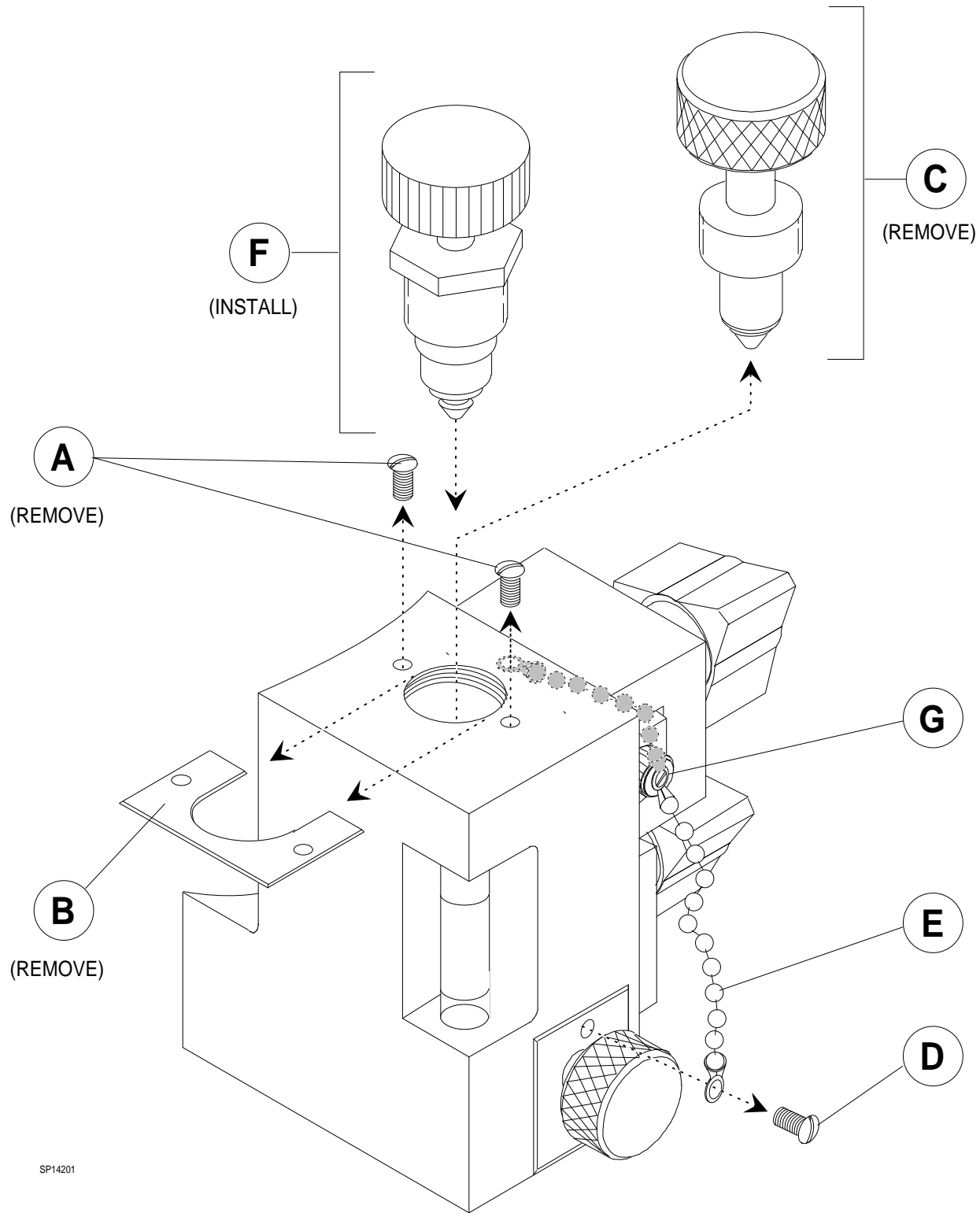

SP14201

Figure 1: VAPORIZER FILL VALVE UPGRADE



NORTH
AMERICAN
DRÄGER

Quality Service for Life®

Technical Service Department
3122 Commerce Drive
Telford, PA 18969
(215) 721-5402
(800) 543-5047
(215) 723-5935 Fax

Part Number: SP00142

Rev: -
Date: June 29, 1994
© N.A.D. Inc. 1994